

SMART TRACK

WEB-BASED APRON TRACKING SYSTEM

Manage, track, and document the inspection of radiation protection garments quickly, easily and completely free.

No hardware or software required.
No registration fees.
No annual fees.
No more hassle.



SmartTrack.InfabCorp.com
805.987.5255



SMART TRACK - *QUICK AND EASY*

We are pleased to introduce our web-based tool that will allow our customers to manage, track, and document the inspections of all of their Radiation Protection Garments.

Smart Track generates reports based on a user's garment inspection history. Users will utilize these reports to give details on inspections during audits from authoritative bodies while also having a much better understanding of their current inventory status.

By inputting your protective garment inventory into Smart Track, you will instantly have an extremely user-friendly way to organize and streamline your inspection process.

Works on PC, Mac and all systems.



FEATURES

- Generates entire inventory, missing, and damaged reports with full inspection histories.
- Users have the ability to filter searches by apron ID, department, and the person assigned.
- Easily request quotes to replace aprons that are damaged or missing.
- Optional utilization of 2D bar-code technology to be able to view and/or inspect specific aprons instantly.

NEW FOR 2015

- **Differentiate between fluoroscopic and visual inspection.**
- **Administrative style accounts**
- **Mobile-friendly: viewable on all computers, tablets, and phones**

BENEFITS

- Streamlines your inspection process.
- Generates simple yet robust reports.
- Digitizes the information that can be accessed from any computer with an Internet connection.
- User-friendly design allows for simple integration.
- Bulk-upload assistance of current inventory is available.



SECURELY WEB-BASED

Smart Track requires no software.
Simply login and register to get started.



ANY APRONS

Track aprons from any manufacturer.

HOW SMART TRACK WORKS

SMART TRACK BASIC INSTRUCTIONS

1 Create Free Account

- Go to SmartTrack.InfabCorp.com
- Click [Sign Up Now](#)
- Enter information and receive confirmation e-mail
- Verify e-mail address & password (case sensitive)

Change Password

- Go to Edit Profile
- Scroll to "New Password"

2 Add Aprons

Select Add Garment

Enter apron information, including:

- Select Inspector (optional admin account)
- I.D. - Unique identifier (mandatory)
- Department
- Assigned to
- Manufacturer (mandatory)
- Garment type (mandatory)
- Core material
- Color/trim
- Monogram
- Manufacturer date

Import Multiple Aprons

- Hover over "Add Garments", select "Import"
- Download the template on the "Import Multiple Aprons" page, fill it out with your current information, and upload it.
- Make any adjustments needed.

Import Garments

Step 1 - Upload File

Maximum file size: 1 MB

Required File Type: Comma Separated Value (.csv)

[Download Template](#) - This is an MS Excel spreadsheet with predefined drop downs for simpler matching. Be sure to save your edited version as a comma separated vale (.csv) before uploading.

Click or drag and drop file here to upload.

[Upload File](#)

3 View Inventory

- Select "Inspect Garments" or "Reports"
- View all aprons on one page or search for a specific apron by field.

Status	I.D.	Name	Inspector	Size	Color	Actions
Inspected	0008	Cap	Dr. O'Connell	Small	Blue/Red	Details
Inspected	10013	APR - Nylon	Dr. O'Connell	Full Wrap	Blue and Green	Details
Inspected	000007			Frontal		Details
Inspected	Sample 1234	Cap-Like	Dr. L.	Small	Blue	Details
Inspected	00003	Cap-Like	Dr. Jane O'Connell	Frontal	Red and Green with a small red and blue logo on the front edge.	Details
Inspected	01024	Ballcap Department	Dr. John Doe	Med and Small	Blue with red trim	Details
Inspected	05676	Cap-Like			Frontal	Details
Inspected	010121	Ballcap Department	Dr. John Doe	Med and Small	Blue with red trim	Details
Inspected	0103	Ballcap Department	Dr. John Doe	Med and Small	Blue with red trim	Details



HOW IT WORKS

4 DOCUMENT INSPECTION RESULTS

1. Click on Inspect to add inspection results.
2. Choose inspection result
 - » Passed Inspection
 - » Failed Inspection
 - » Missing
3. Complete Inspection

ARCHIVE APRONS

Select "Archive" to indicate apron is no longer in use.

Inspect Garments [Add New Garment](#)

Select to view garments for: **All Inspectors** [assign garments to inspectors](#)

Type of inspection today: Fluoroscopic Visual

[Save All](#) [Scan](#)

[expand all](#) [clear filters](#)

Passed	Failed	Missing	Next Inspection	Status	Apron Id	Department	Type	Color	
<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	in 6 months	Passed	S0008	Cath	Breast Shield	Blue/Red	details

Inspection Notes:

Inspection Type: Fluoroscopic Visual

Inspector: Jones, Tom5
Assigned: Dr. O'Connell
Manufacturer: AADCO
Size: Custom
Core: KIARMOR
Protection: 0.25 mm
Monogram: Dr. O'CI
Purchased/Manufactured: Feb 2014
Last Updated: 10 Dec 2014 12:12 PM

Last Inspected: 10 Dec 2014 12:11 PM
Last Inspection Type: Fluoroscopic
Last Inspected By: User, Demo
Last Inspection Notes:

[edit](#) | [archive](#) | [report](#) | [delete](#)

5 VIEW REPORTS

Select Reports or rollover Reports to see additional report options: Full Apron Inventory, Damaged Aprons, Missing Aprons and Aprons Removed from Service.

All Garments

Select to view garments for: **All Inspectors** [assign garments to inspectors](#) [Scan](#)

[expand all](#) [clear filters](#)

Apron Id	Department	Type	Last Inspection	Status	
<input type="text" value="Search..."/>	<input type="text" value="Search..."/>	<input type="text" value="Search..."/>	<input type="text" value="Search..."/>	<input type="text" value="Search..."/>	
ID2345	Cath Lab	Half Apron	24 Nov 2014 9:05 AM	Archived	view report
S0008	Cath	Breast Shield	10 Dec 2014 12:11 PM	Passed	view report
test 123	RSNA	Full Wrap	10 Dec 2014 12:11 PM	Passed	view report
S000807		Frontal	10 Dec 2014 12:11 PM	Failed	view report
sample 1234	Cath Lab	Skirt	10 Dec 2014 12:11 PM	Passed	view report

Single Garment

ID ID2345
 Department: Cath Lab
 Assigned To: Dr. Jane Doe
 Manufacturer: Infab
 Garment Type: Half Apron
 Core Material: Lead Free
 Size: Custom
 Protection Level: 1.0 mm
 Color/Trim: All black
 Monogram: Dr. Jane
 Manufactured Date: Dec 2013
[edit](#)

Inspection History

Inspection Date	Status	Inspected By	Inspection Type	Notes
24 Nov 2014 9:05 AM	Archived	User, Demo		
02 Dec 2014 4:20 PM	Passed			Garment imported: 2014-11-21



SMART TRACK


FREE MOBILE APP

Simply download the app for free, login to your Smart Track account and scan and track all of your aprons quickly, easily and completely free.

The Smart Track mobile app makes it easy to scan and track all of the aprons in a department or facility. The Smart Track app has all the features of the web-based version of Smart Track **but also includes a built in scanner. The Smart Track app requires no sign up fee and no annual fee making it the easiest and best way to track radiation protection garments.**

INSTALL APP

1. Search the iTunes app store for "Infab Smart Track"
2. Download free app seen below
3. Login and start tracking

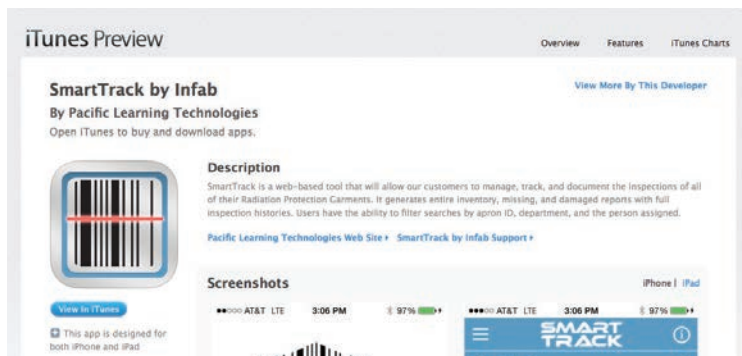



SMART TRACK RETROFIT TAG KIT

ST-OPT-TAG

Smart Track Scanner and Smart Track Tags with Backing

Includes scheduled training webinars




1 SMART TRACK Scanner



200 Annual Inspection Stickers



200 SMART TRACK 2D bar code tags with permanent security-style backing

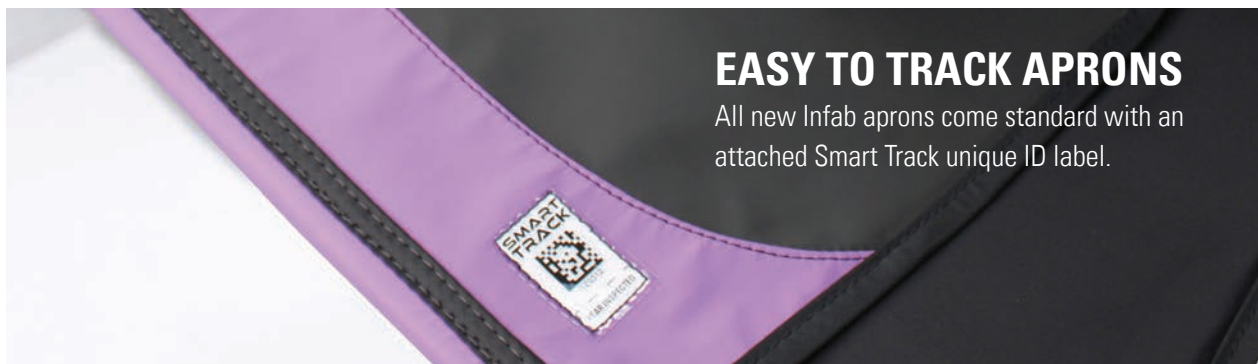

EASY TO TRACK APRONS

All new Infab aprons come standard with an attached Smart Track unique ID label.

CUSTOMER SATISFACTION



“ *Smart Track makes it simple to manage my suite of lead aprons. We track over 100 aprons - it saves time and headaches.* ”

Rich Norman
Lead Tech, Los Robles Hospital
Thousand Oaks, CA

“ *The SMART TRACK Apron Inventory system has revolutionized the way lead apron quality control testing and inventory management is performed at our facility. It is a superior, comprehensive, and user-friendly product. Individual bar codes make tracking aprons quick and easy. It is the best system I have had the opportunity to work with. —* ”

Doug Beardsall MRT (R) CTIC
Senior Technologist QC, Safety, and
Education, Royal Victoria Regional
Health Center-- Barrie, Ontario

“ *Tracking lead aprons and making sure that they have been inspected is a challenge for any large radiology facility. Of the many apron-tracking systems I reviewed, I knew that the “SMART TRACK” Apron inventory system would work for our facility and meet new regulatory standards set by “The Joint Commission”... The system is user friendly. All I know is that anyone who has the burden of tracking aprons needs to get “Smart” about tracking their apron inventory”.* ”

Roseann Zeldin, RT, (R) (M), CNMT
Assistant Radiation Safety Officer,
Department of Radiology & Imaging,
Hospital for Special Surgery-- New York
City, New York



SmartTrack.InfabCorp.com



Contact Us
SmartTrack@Infabcorp.com
M-F 6:30am-5pm PDT
info@InfabCorp.com
www.InfabCorp.com

Infab Corporation
1040 Avenida Acaso
Camarillo, CA, USA 93012
Phone: 805-987-5255
Fax: 805-482-8424